


Fibre photometry 1.07 installation notes

1. Uninstall Fibre photometry 1.06 using Windows system settings. Keep NIDAQ

2. Delete all settings files for previous software version. These can be found wherever the raw data is exported or has been previously exported to. Best to search for “fibre_photometry.settings” and delete all. See screenshots on left.

Restart computer

3. Install Fibre photometry 1.07
Do not install NIDAQMX again

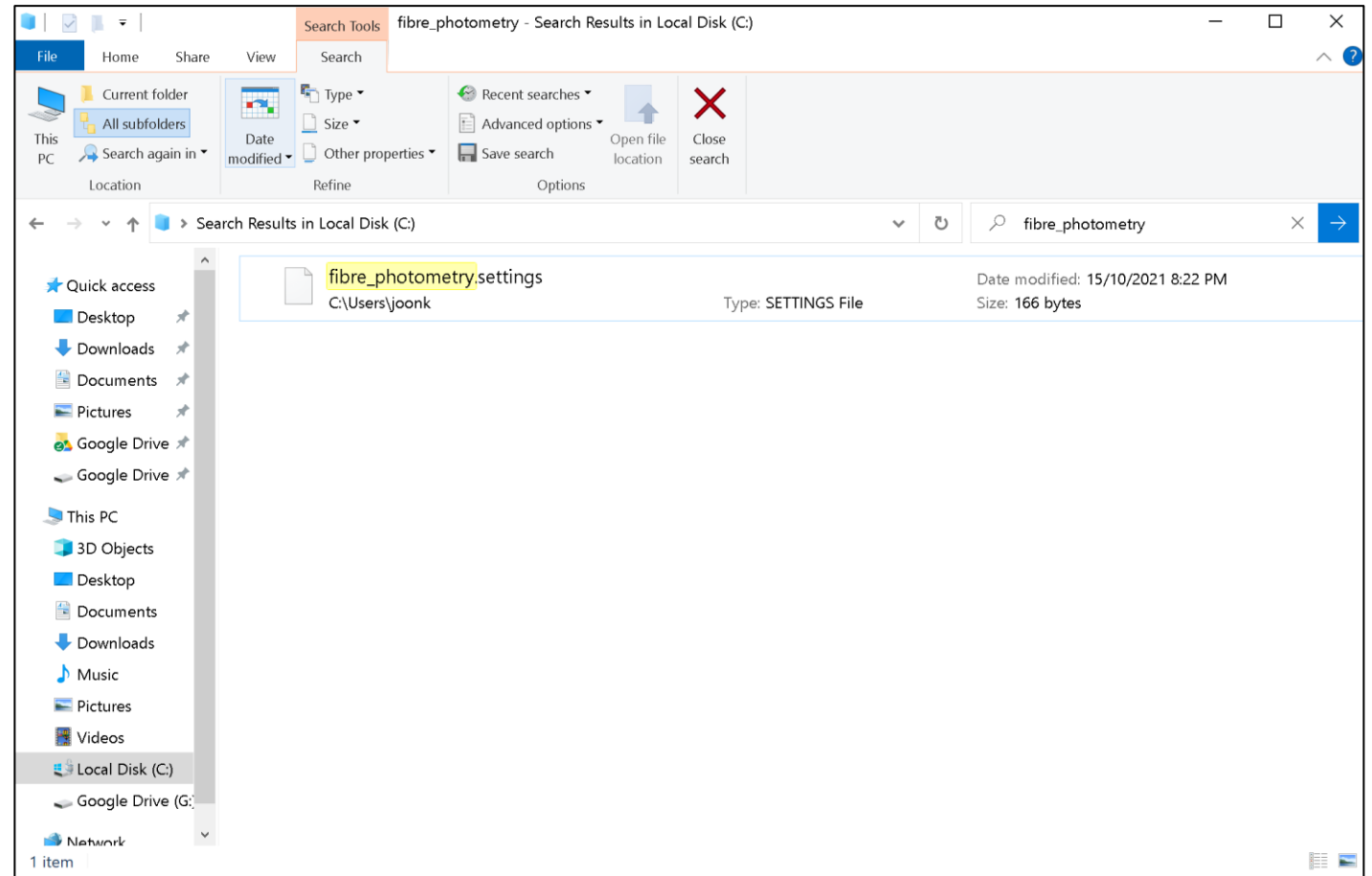

-1.14-windows	22/08/2021 8:09 PM	XML Document	7 KB
fibre_photometry.settings	15/10/2021 8:22 PM	SETTINGS File	1 KB
NTUSER.DAT	12/10/2021 9:50 PM	DAT File	8,192 KB